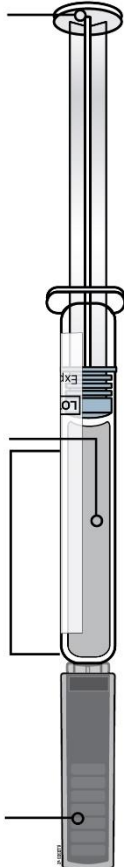
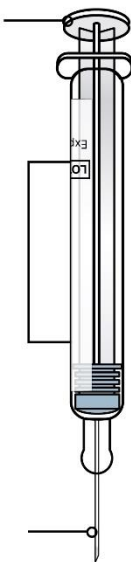



Instructions for use:
Repatha® single use pre-filled syringe

Guide to parts	
Before use	After use
<p>Plunger rod</p>  <p>Medicine</p> <p>Syringe barrel</p> <p>Grey needle cap on</p>	<p>Used plunger</p>  <p>Used syringe barrel</p> <p>Used Needle</p> <p>Grey needle cap off</p>
	Needle is inside.

Important

Before you use a single use Repatha pre-filled syringe, read this important information:

- **Do not** freeze or use the Repatha pre-filled syringe if it has been frozen.
- **Do not** use the Repatha pre-filled syringe if the packaging is open or damaged.
- **Do not** use the Repatha pre-filled syringe if it has been dropped onto a hard surface. Part of the syringe may be broken even if you cannot see the break. Use a new Repatha pre-filled syringe.
- **Do not** remove the grey needle cap from the Repatha pre-filled syringe until you are ready to inject.

Step 1: Prepare

A	Remove the Repatha pre-filled syringe carton from the refrigerator and wait 30 minutes.
----------	--

Wait at least 30 minutes for the pre-filled syringe in the carton to naturally reach room temperature before injecting.

Check that the name Repatha appears on the carton label.

- **Do not** try to warm the Repatha pre-filled syringe by using a heat source such as hot water or microwave.
- **Do not** leave the Repatha pre-filled syringe exposed to direct sunlight.
- **Do not** shake the Repatha pre-filled syringe.


B	Gather all materials needed for your injection.
----------	--

Wash your hands thoroughly with soap and water.

On a clean, well-lit, flat work surface, place:

- One Repatha pre-filled syringe in its tray.
- Alcohol wipes.
- Cotton ball or gauze pad.
- Plaster.
- Sharps disposal container.
- **Do not** use if expiry date on the Repatha pre-filled syringe carton has passed.

C	Choose your injection site.
----------	------------------------------------

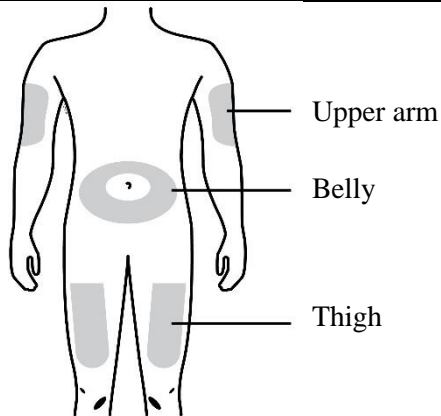


You can use:

- Thigh.
- Belly, except for the 2 inches (5 centimetres) around the belly button.
- Outer area of upper arm (only if someone else is giving you the injections).
- **Do not** choose an area where the skin is tender, bruised, red, or hard. Avoid injecting into areas with scars or stretch marks.



Choose a different site each time you give yourself an injection. If you need to use the same injection site, just make sure it is not the same spot on that site you used last time.

D	Clean your injection site.
----------	-----------------------------------

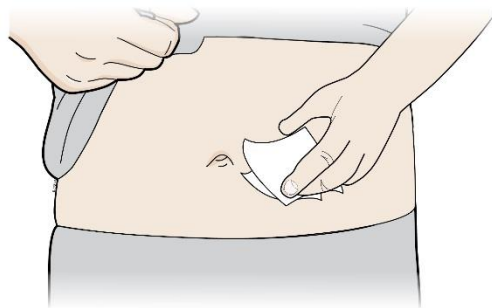




Clean your injection site with an alcohol wipe. Let your skin dry before injecting.

- **Do not** touch this area of skin again before injecting.

E	Remove the pre-filled syringe from the tray.
----------	---

Turn tray over

Gently Press

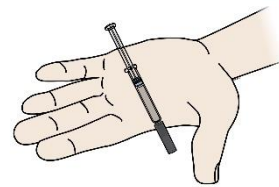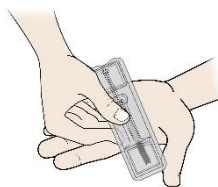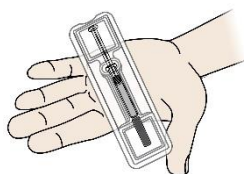

To remove:

- Peel the paper off of the tray.
- Place the tray on your hand.
- Turn the tray over and gently press the middle of the tray's back to release the syringe into your palm.
- If the pre-filled syringe does not release from the tray, gently press on the back of tray.
- **Do not** pick up or pull the pre-filled syringe by the plunger rod or grey needle cap. This could damage the syringe.
- **Do not** remove the grey needle cap from the pre-filled syringe until you are ready to inject.



Always hold the pre-filled syringe by the syringe barrel.

F	Inspect medicine and syringe.
----------	--------------------------------------

Plunger rod

Syringe barrel

Syringe label with expiration date

Grey needle cap on



Medicine



Always hold the pre-filled syringe by the syringe barrel.

Check that:


- The name Repatha appears on the pre-filled syringe label.
- The medicine in the pre-filled syringe is clear and colourless to slightly yellow.
- **Do not** use the pre-filled syringe if any part of the pre-filled syringe appears cracked or broken.
- **Do not** use the pre-filled syringe if the grey needle cap is missing or not securely attached.
- **Do not** use the pre-filled syringe if the medicine is discoloured or contains large lumps, flakes or coloured particles.
- **Do not** use the pre-filled syringe if the expiration date on the pre-filled syringe has passed.

Step 2: Get ready	
A	Carefully pull the grey needle cap straight out and away from your body. Do not leave the grey needle cap off for more than 5 minutes. This can dry out the medicine.
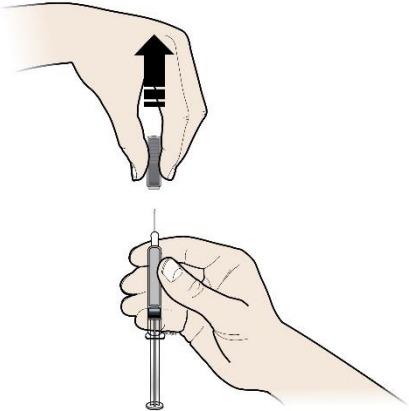
<p>1.</p>  <p>It is normal to see a drop of medicine at the end of the needle.</p>	<p>2.</p>  <p>Immediately place the cap in the sharps disposal container.</p>
<ul style="list-style-type: none"> • Do not twist or bend the grey needle cap. This can damage the needle. • Do not put the grey needle cap back onto the pre-filled syringe. 	

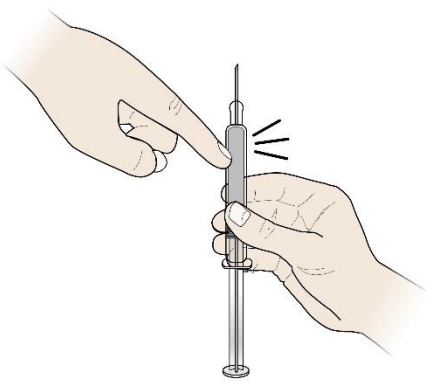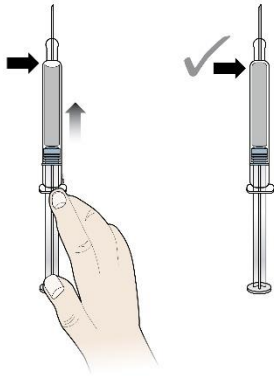B	Remove the air bubble/gap.
<p>You may notice an air bubble/gap in the Repatha pre-filled syringe.</p> <p>If you notice an air bubble/gap:</p> <ul style="list-style-type: none"> • Hold the pre-filled syringe with the needle facing up. • Gently tap the syringe barrel with your fingers until the air bubble/gap rises to the top of the syringe. • Slowly and gently push the plunger rod up to get the air out of the pre-filled syringe. Be very careful not to push out any medicine. <div style="display: flex; justify-content: space-around; align-items: center;">   </div> <ul style="list-style-type: none"> • Do not tap the syringe needle. 	

C **PINCH** your injection site to create a firm surface.

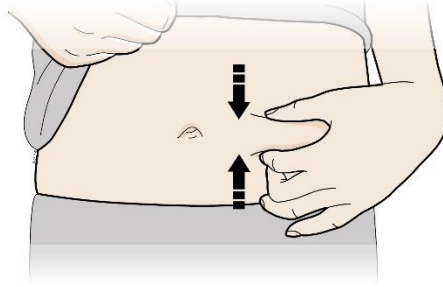


Pinch skin firmly between your thumb and fingers, creating an area about 2 inches (5 centimetres) wide.



It is important to keep the skin pinched while injecting.


Step 3: Inject


A **Hold the PINCH.** Insert the needle into the skin using a 45 to 90 degree angle.

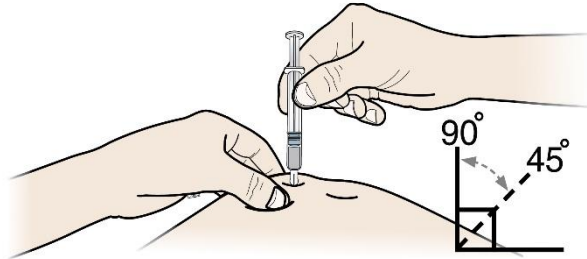


- **Do not** place your finger on the plunger rod while inserting the needle.

B Using slow and constant pressure, **PUSH** the plunger rod all the way down until the syringe is empty.

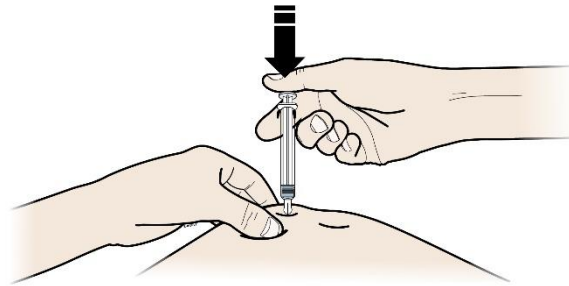

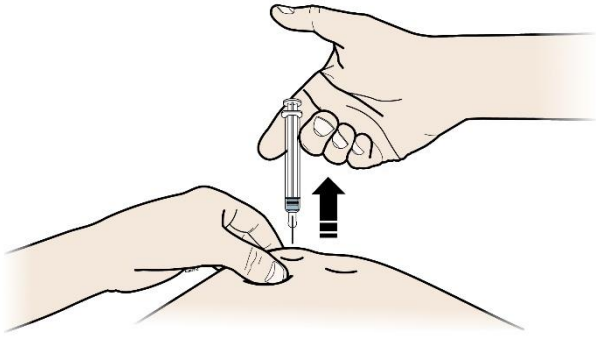C	When done, RELEASE your thumb, and gently lift the syringe off skin.
	
<ul style="list-style-type: none"> • Do not put the grey needle cap back onto the used syringe. 	

Step 4: Finish	
A	Immediately place the used syringe in a sharps disposal container.
	
<p>Talk with your healthcare provider about proper disposal. There may be local guidelines for disposal.</p> <ul style="list-style-type: none"> • Do not reuse the used syringe. • Do not use any medicine that is left in the used syringe. • Do not recycle the syringe or the sharps disposal container or throw it into household rubbish. <p>! Keep the used syringe and sharps container out of the sight and reach of children.</p>	

B	Examine the injection site.
<p>If there is blood, press a cotton ball or gauze pad on your injection site. Apply a plaster if needed.</p> <ul style="list-style-type: none"> • Do not rub the injection site. 	

January 2021
SGREPPFSIFU03

AMGEN®

Repatha® is a registered trademark owned or licensed by Amgen Inc., its subsidiaries, or affiliates.