


Instructions for Use

Emgality Solution for Injection in Pre-filled Syringe 100 mg/ml

Galcanezumab

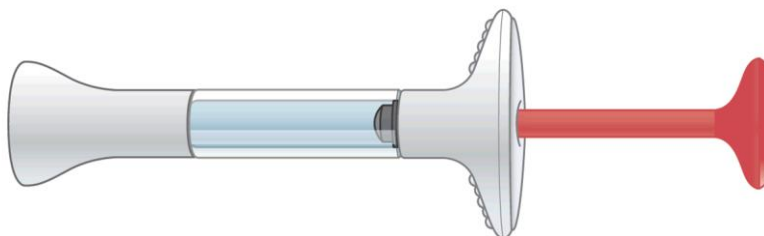

For subcutaneous injection only.

Before you use the Emgality pre-filled syringe, read and carefully follow all the step-by-step instructions.


Important Information

- Your healthcare provider or nurse should show you how to prepare and inject Emgality using the pre-filled syringe. Do not inject yourself or someone else until you have been shown how to inject Emgality.
- Keep this Instructions for Use and refer to it as needed.
- Each Emgality pre-filled syringe is for **one-time use only**. **Do not** share or reuse your Emgality pre-filled syringe. You may give or get an infection.
- Your healthcare provider may help you decide where on your body to inject your dose. You can also read the “**Choose your injection site**” section of these instructions to help you choose which area can work best for you.
- If you have vision or hearing problems, **do not** use Emgality pre-filled syringe without help from a caregiver.
- See “**Storage and Handling Information**” for important storage information.

INSTRUCTIONS FOR USE

Before you use the Emgality pre-filled syringe, read and carefully follow all the step-by-step instructions.

Parts of the Emgality Pre-filled Syringe

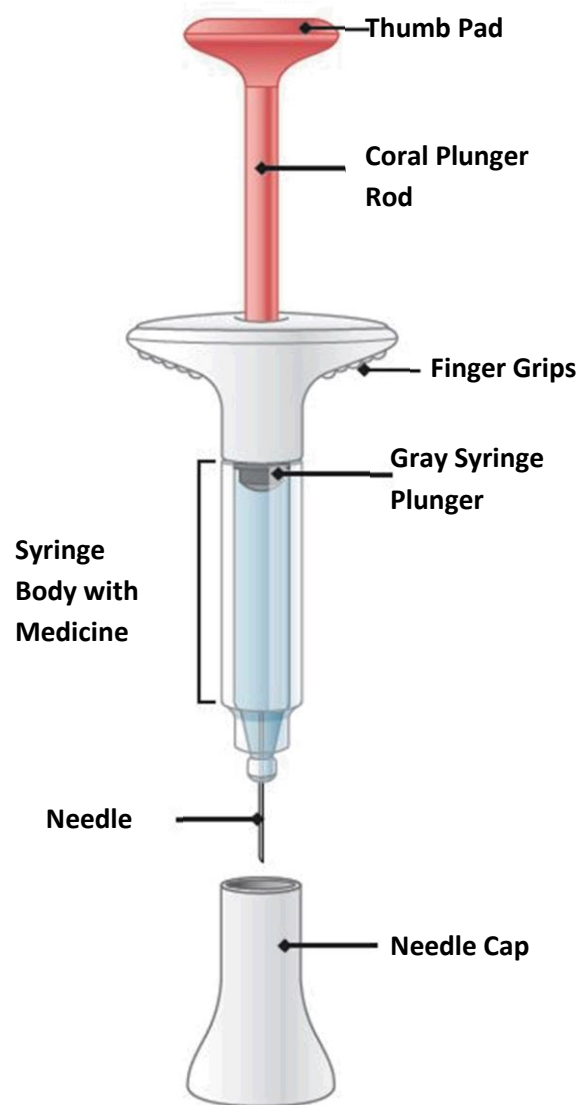

Before You Get Started

Take the Pre-filled Syringe from the refrigerator

Take 3 Emgality pre-filled syringes from the refrigerator.

Check your prescription.

- Emgality comes as a single-dose pre-filled syringe.
- You will need 3 pre-filled syringes for your each dose.

Leave the needle caps on until you are ready to inject.

Leave the pre-filled syringes at room temperature for 30 minutes before injecting.

Do not microwave the pre-filled syringe, run hot water over it, or leave it in direct sunlight.

Do not shake.

Gather Supplies

For each injection you will need:

- 1 alcohol wipe
- 1 cotton ball or piece of gauze
- 1 sharps disposal container. See “After You Inject Your Medicine.”


Inspect the Pre-filled Syringe and the medicine

Make sure you have the right medicine. The medicine inside should be clear. Its color may be colorless to slightly yellow.

Do not use the pre-filled syringe, and throw away (dispose of) as directed by your healthcare provider or pharmacist if:

- it looks damaged
- the medicine is cloudy, is discolored, or has small particles
- the Expiration Date (Exp.) printed on the label has passed
- the medicine is frozen


Expiration Date

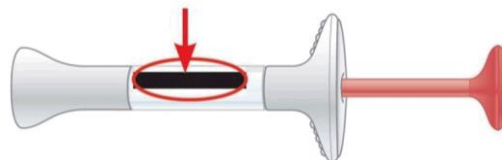

Prepare for injection

Wash your hands with soap and water before you inject your Emgality. Make sure a sharps disposal container is close by.

Choose your injection site

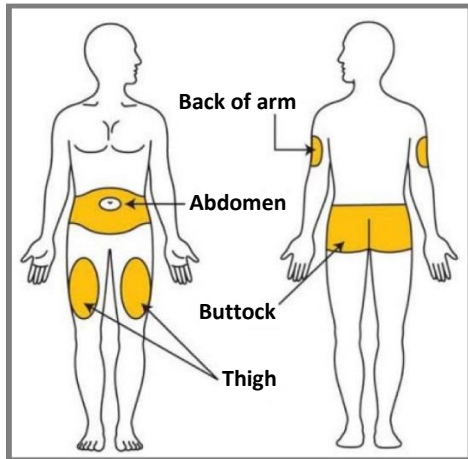


Your healthcare provider can help you choose the injection site that is best for you.

- **You** may inject the medicine into your stomach area (abdomen). Do not inject within 2 inches of the belly button (navel).
- **You** may inject the medicine into the front of your thighs. This area should be at least 2 inches above the knee and 2 inches below the groin.
- **Another person** may give you the injection in the back of your upper arm or buttocks.
- **Do not** inject in the exact same spot. For example, if your first injection was in your abdomen, your next injection could be in another area of your abdomen.
- **Do not** inject into areas where the skin is tender, bruised, red, or hard.
- **Clean your injection site with an alcohol wipe. Let the injection site dry before you inject.**


1 Uncap

- **Leave the needle cap on until you are ready to inject.**
- Pull the needle cap off and throw it away in your household trash.
- **Do not** put the needle cap back on – you could damage the needle or stick yourself by accident.
- **Do not** touch the needle.

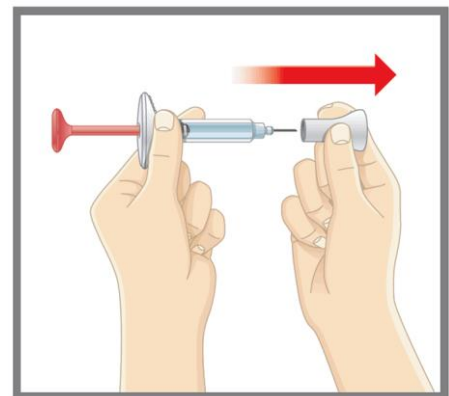


2 Insert

- Gently pinch and hold a fold of skin where you will inject.
- Insert the needle at a 45-degree angle.

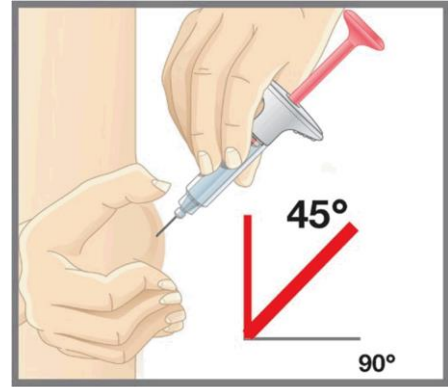


3 Inject

- Slowly push on the thumb pad to push the plunger all the way in until all the medicine is injected.
- The gray syringe plunger should be pushed all the way to the needle end of the syringe.

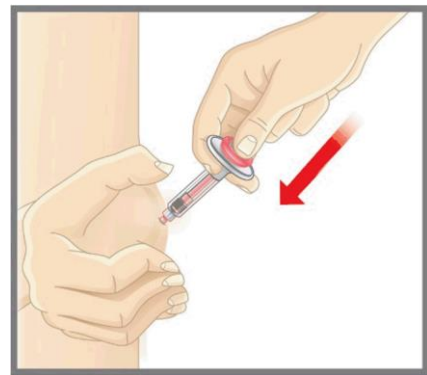

- You should see the coral plunger rod show through the syringe body when the injection is complete as shown.
- Remove the needle from your skin and gently let go of your skin.
- If you have bleeding at the injection site, press a cotton ball or gauze over the injection site. Do not rub the injection site.
- **Do not** put the needle cap back on the pre-filled syringe.

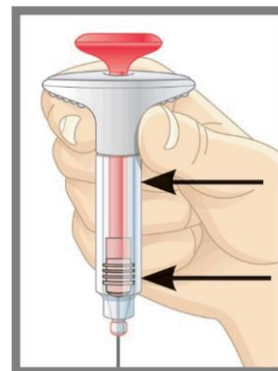

Coral plunger rod

Gray syringe plunger

After You Inject Your Medicine

Throw away the used pre-filled syringe

- Put the used Emgality pre-filled syringe in a sharps disposal container right away after use.
Do not throw away (dispose of) the Emgality pre-filled syringe in your household trash.



- If you do not have a sharps disposal container, ask your doctor, pharmacist or nurse where you can get one.
- When your sharps disposal container is almost full, dispose as directed by your doctor, pharmacist or nurse.
- **Do not** recycle your used sharps disposal container.

For each of the 3 injections, repeat all instructions with a new pre-filled syringe.

Commonly Asked Questions

Q. What if I see air bubbles in my Emgality pre-filled syringe?

A. It is normal to have air bubbles in the pre-filled syringe. Emgality is injected under your skin (subcutaneous injection), so these air bubbles will not harm you.

Q. What if there is a drop of liquid on the tip of the needle when I remove the needle cap?

A. It is okay to see a drop of liquid on the tip of the needle.

Q. What if I cannot push in the plunger?

A. If the plunger is stuck or damaged:

- **Do not** continue to use the syringe
- Remove the needle from your skin
- Dispose of the syringe and get a new one

Q. What if there is a drop of liquid or blood on my skin after my injection?

A. This is normal. Press a cotton ball or gauze over the injection site. Do not rub the injection site.

Q. How can I tell if my injection is complete?

A. When your injection is complete:

- The coral plunger rod should show through the body of the syringe.
- The gray syringe plunger should be pushed all the way to the needle end of the syringe.

If you have more questions about how to use the Emgality pre-filled syringe, contact your doctor, pharmacist or nurse.

Storage and Handling Information

- Store your pre-filled syringes in the refrigerator between 2°C to 8°C.
- Your pre-filled syringes may be stored outside of the refrigerator in the original carton at temperatures up to 30°C for up to 7 days. After storing out of the refrigerator, **do not** place Emgality back in the refrigerator. **Do not** store above 30°C.
- **Do not** freeze your pre-filled syringes.
- Keep your pre-filled syringes in the carton it comes in to protect your pre-filled syringes from light until time of use.
- **Do not** shake your pre-filled syringes.
- Throw away your pre-filled syringes if any of the above conditions are not followed.
- **Keep your pre-filled syringes and all medicines out of the reach of children.**

Read the full Prescribing Information for Emgality inside this box to learn more about your medicine.

Product Owner

Eli Lilly and Company, Indianapolis, IN 46285, USA

Date of revision of text: 18 Jan 2021